



Nilocin Wart & Verruca Treatment Gel

UK



For external use only. Keep out of sight and reach of children



Nilocin® contains a corrosive agent.



Read this instruction carefully prior to use

INDICATION

Nilocin® Wart Gel is indicated for the treatment of warts and verrucas found on the feet, hands, body arms and legs. The product can be used by adults and children from 4 years old.



WART AND VERRUCA

Warts and verruca (hereafter defined as warts) are local growths in the skin caused by the human papilloma virus (HPV). They may be found alone or in clusters. There exist different types of warts, most commonly found on hands, feet, elbows and knees. Warts are contagious and spread by contact of the virus with skin.

NILOCIN®

Nilocin Wart treatment is a CE marked Medical Device, for topical administration directly on warts. The pack comprises a glass bottle with 4 ml treatment gel, (enough to treat 10 warts) and this instruction leaflet.

The gel contains a keratolytic compound, which breaks down tissue locally and thereby eliminates warts.



Warnings & Precautions

For external use only. Keep out of reach of children. The product contains a corrosive agent.

- Do NOT use Nilocin® for genital warts or in the facial area.
- Do not use Nilocin® if you are known to be hypersensitive or allergic to the ingredients.
- When applying, avoid contact with healthy skin surrounding the wart. If the gel makes contact with healthy skin, rinse the area with soap and plenty of water.
- Do not expose the treated area to sunlight, until the skin has returned to normal.
- Avoid contact with i.e. jewelry, furniture or clothes, as the gel contains a corrosive agent
- Do not apply Nilocin® to sore or bleeding warts, irritated, inflamed, blistered or damaged skin. Pause until the skin has recovered
- Do not use on warts with hairs, birthmarks, dark moles, pigmented skin or in case you have doubts if the skin condition is a wart.
- Nilocin is not suitable for use if you have diabetes or poor blood circulation.

• it is better to apply too little gel if in doubt. Over application can result in side effects such as wounds, blisters, inflammation, tissue death (necrosis), scars, and chemical burns

• Consult your doctor if the treated area becomes red and painful or gets inflamed.

• Single patient - multiple use

• Do not use in combination with other wart removal therapies.

• Do not peel of the scab, as it needs to fall off by itself

• Do not use for children under 4 years of age, unless directed by doctor.

• If you are pregnant, breastfeeding or have a skin disease, it is recommended to consult your doctor prior to use.

SHOULD THE PRODUCT ACCIDENTALLY COME INTO CONTACT WITH THE EYES OR MUCOSAL MEMBRANES, RINSE WITH PLENTY OF WATER FOR 15 MINUTES AND SEEK MEDICAL ADVICE

HOW TO USE NILOCIN

1 Remove the cap. Push down and turn the cap counterclockwise.

2 Apply gel to the wart with the applicator, twice a day for 4 days. Larger warts might need more than 1 application to cover the entire wart. For children with small warts, treatment schedule is reduced to one application a day

3 Recap the bottle firmly after use

4 Let the gel dry for 10-15 minutes

5 Repeat treatment for 4 days. If the wart has not disappeared completely, repeat the treatment after a 4 days break. This treatment schedule can be repeated up to 4 times.

Make sure to cover the whole wart and take care not to touch or spill onto surrounding skin.

Should any discomfort or pain occur, stop application. If it burns, sting, or if gel is applied to healthy skin, rinse with soap and cold water.

Do not exceed the number of recommended applications.

POSSIBLE SIDE EFFECTS

It is important to use the product as instructed. Do not over apply as this can cause skin damage! (please see PRECAUTIONS above).

A mild tingling or burning sensation can in some cases be felt for 10-15 minutes. If gel is over applied side effects can occur. Pause treatment until the skin has fully recovered.


Serious incidents should be reported to the manufacturer and competent authority

STORAGE

Store between 2 to 25 degree Celsius and out of direct sunlight, enclosed in the original packaging and out of sight and reach of children. Do not use after the expiry date, or 6 months after removing the cap from the bottle for the first time. Carefully recap the bottle after use.

The empty product can be disposed with.

Unigroup ApS, Diplomvej 373, DK-2800 Lyngby, Denmark
Nilocin is a registered trademark of Unigroup ApS.

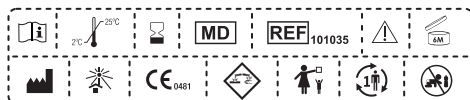

TECHNICAL DATA – Product name: nilocin Product: Wart & Verruca Treatment gel **Ingredients (INCI):** Trichloroacetic acid (TCR 40%), Aqua, Glycerin, Cellulose, Cellulose Gum, Benzoic Acid, Sodium Benzoate **Content:** 4 ml, sufficient to treat minimum 10 warts. **Registration:** Nilocin® is classified and registered as a Medical Device Class IIa. **Storage:** Store at 2-25 degree Celsius, in the original packaging, out of direct sunlight and out of reach and sight of children. Not to be used after the expiry date. **Manufacturer:** Unigroup PPS, Diplomvej 373, DK-2800 Lyngby, Denmark. www.nilocin.com. Information updated: 2020-02

UG201908-101035-UK



Nilocin Wart & Verruca Treatment Gel

UK



For external use only. Keep out of sight and reach of children



Nilocin® contains a corrosive agent.



Read this instruction carefully prior to use

INDICATION

Nilocin® Wart Gel is indicated for the treatment of warts and verrucas found on the feet, hands, body arms and legs. The product can be used by adults and children from 4 years old.



WART AND VERRUCA

Warts and verruca (hereafter defined as warts) are local growths in the skin caused by the human papilloma virus (HPV). They may be found alone or in clusters. There exist different types of warts, most commonly found on hands, feet, elbows and knees. Warts are contagious and spread by contact of the virus with skin.

NILOCIN®

Nilocin Wart treatment is a CE marked Medical Device, for topical administration directly on warts. The pack comprises a glass bottle with 4 ml treatment gel, (enough to treat 10 warts) and this instruction leaflet.

The gel contains a keratolytic compound, which breaks down tissue locally and thereby eliminates warts.



Warnings & Precautions

For external use only. Keep out of reach of children. The product contains a corrosive agent.

- Do NOT use Nilocin® for genital warts or in the facial area.
- Do not use Nilocin® if you are known to be hypersensitive or allergic to the ingredients.
- When applying, avoid contact with healthy skin surrounding the wart. If the gel makes contact with healthy skin, rinse the area with soap and plenty of water.
- Do not expose the treated area to sunlight, until the skin has returned to normal.
- Avoid contact with i.e. jewelry, furniture or clothes, as the gel contains a corrosive agent
- Do not apply Nilocin® to sore or bleeding warts, irritated, inflamed, blistered or damaged skin. Pause until the skin has recovered
- Do not use on warts with hairs, birthmarks, dark moles, pigmented skin or in case you have doubts if the skin condition is a wart.
- Nilocin is not suitable for use if you have diabetes or poor blood circulation.

• it is better to apply too little gel if in doubt. Over application can result in side effects such as wounds, blisters, inflammation, tissue death (necrosis), scars, and chemical burns

• Consult your doctor if the treated area becomes red and painful or gets inflamed.

• Single patient - multiple use

• Do not use in combination with other wart removal therapies.

• Do not peel of the scab, as it needs to fall off by itself

• Do not use for children under 4 years of age, unless directed by doctor.

• If you are pregnant, breastfeeding or have a skin disease, it is recommended to consult your doctor prior to use.

SHOULD THE PRODUCT ACCIDENTALLY COME INTO CONTACT WITH THE EYES OR MUCOSAL MEMBRANES, RINSE WITH PLENTY OF WATER FOR 15 MINUTES AND SEEK MEDICAL ADVICE

HOW TO USE NILOCIN

1 Remove the cap. Push down and turn the cap counterclockwise.

2 Apply gel to the wart with the applicator, twice a day for 4 days. Larger warts might need more than 1 application to cover the entire wart. For children with small warts, treatment schedule is reduced to one application a day

3 Recap the bottle firmly after use

4 Let the gel dry for 10-15 minutes

5 Repeat treatment for 4 days. If the wart has not disappeared completely, repeat the treatment after a 4 days break. This treatment schedule can be repeated up to 4 times.

Make sure to cover the whole wart and take care not to touch or spill onto surrounding skin.

Should any discomfort or pain occur, stop application. If it burns, sting, or if gel is applied to healthy skin, rinse with soap and cold water.

Do not exceed the number of recommended applications.

POSSIBLE SIDE EFFECTS

It is important to use the product as instructed. Do not over apply as this can cause skin damage! (please see PRECAUTIONS above).

A mild tingling or burning sensation can in some cases be felt for 10-15 minutes. If gel is over applied side effects can occur. Pause treatment until the skin has fully recovered.


Serious incidents should be reported to the manufacturer and competent authority

STORAGE

Store between 2 to 25 degree Celsius and out of direct sunlight, enclosed in the original packaging and out of sight and reach of children. Do not use after the expiry date, or 6 months after removing the cap from the bottle for the first time. Carefully recap the bottle after use.

The empty product can be disposed with.

Unigroup ApS, Diplomvej 373, DK-2800 Lyngby, Denmark
Nilocin is a registered trademark of Unigroup ApS.

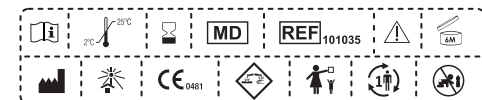

TECHNICAL DATA – Product name: nilocin Product: Wart & Verruca Treatment gel **Ingredients (INCI):** Trichloroacetic acid (TCR 40%), Aqua, Glycerin, Cellulose, Cellulose Gum, Benzoic Acid, Sodium Benzoate **Content:** 4 ml, sufficient to treat minimum 10 warts. **Registration:** Nilocin® is classified and registered as a Medical Device Class IIa. **Storage:** Store at 2-25 degree Celsius, in the original packaging, out of direct sunlight and out of reach and sight of children. Not to be used after the expiry date. **Manufacturer:** Unigroup PPS, Diplomvej 373, DK-2800 Lyngby, Denmark. www.nilocin.com. Information updated: 2020-02

UG201908-101035-UK